Please read the instructions before using your Panthera D-SAD Classic or D-SAD X3.

Introduction

The D-SAD Classic and D-SAD X3 are devices for the treatment of snoring and mild to moderate obstructive sleep apnea (OSA) in adults. The device requires a prescription from a physician. The Panthera D-SAD devices are made to your precise measurements and are single patient use. How do these devices work? When worn during sleep, the devices are positioned on the teeth to move the lower jaw forward. By advancing the jaw, the tongue is also moved forward to help keep the upper airway open. Opening the airway eases the airflow, reducing snoring and OSA.

Indications for Use

The Panthera D-SAD devices are intended to reduce or alleviate snoring and mild to moderate obstructive sleep apnea (OSA) in adults.

Important to remember

- Always wear both upper and lower splints together. Before each use, check the device for:
- Possible damage: the device should not be used if there is visible wear or damage.
- Cleanliness: if the device is not clean, it should be cleaned as per the cleaning instructions.

Contraindications

- The Panthera devices are contraindicated for patients who:
- Are under 18 years of age.
 - Have loose teeth.
 - Advanced periodontal disease.
 - Severe respiratory disorders other than OSA.
 - Central sleep apnea.

Precautions

- The dental sleep specialist should consider the medical history of the patients, including history of asthma, breathing or respiratory disorders, or other relevant health problems, and refer the patient to the appropriate healthcare provider before prescribing the device.
- The Panthera devices must be used during sleep at night or when taking a nap. Do not eat while wearing the device. It is possible to drink water.
- It is recommended to consult a dental professional for an annual checkup.

Possible Side Effects


- The use of the Panthera devices may cause:
- Tooth movement or changes in dental occlusion.
 - Gingival or dental soreness.
 - Pain or soreness of the temporomandibular joint.
 - Obstruction of oral breathing.
 - Excessive salivation or mouth dryness.
 - Any abnormal incident related to wearing the device must be reported to Panthera Dental and for EU patients, must also be reported to the competent authority in your country of residence.
 - Should you experience significant side effects, contact your dental practitioner.

Package Contents

- D-SAD Classic with rod kit OR D-SAD X3 with clip kit.
- Cleaning bath.
- Storage container.
- Instructions for use.
- Optional: Orthodontic elastic bands.
- Optional: Patient dental model.

Parts

- Upper splint and lower splint with fixed pivot point (anchor).
- A rod kit for jaw advancement. The rods come in different lengths; the dentist will determine the best length for you.

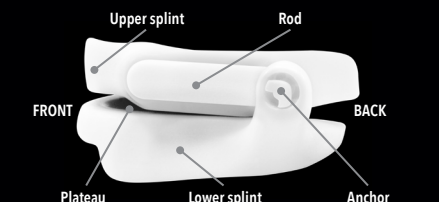
Note: The design of the device may vary depending on the specifics requested by your dental professional.

Rod Attachment
Rods connect the upper and lower splints (hinge) and allow to move the lower jaw forward.



- Attach one rod on each side of the lower splint into the anchor cutout, the fixed anchor point will match the cutout on the rod.
- Press firmly to secure. The rod number should be facing inward and the unattached end should be directed to the back of the splint.

Attach the upper and lower splints on one side first.

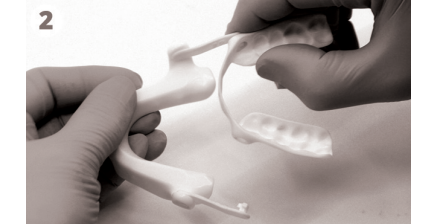- | | |
|---|--|
| Prescription only <ul style="list-style-type: none">Sous prescription seulementVerschreibungspflichtigSolo con ricetta medicaSolo su prescrizioneUitsluitend op recept | Authorized European representative <ul style="list-style-type: none">Représentant autorisé européenBevollmächtigter europäischer VertreterRepresentante europeo autorizadoRappresentante europeo autorizzatoBevoegde Europees vertegenwoordiger |
|---|--|

- To attach the upper splint, position it with the "teeth" side showing (facing upward); attach one side of the upper splint to the lower splint with the rod inserted into the cutout.

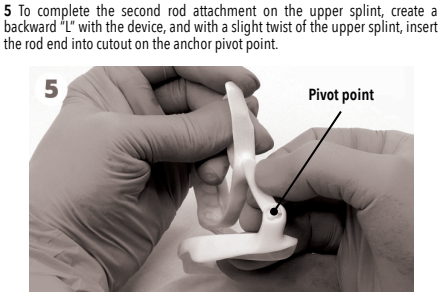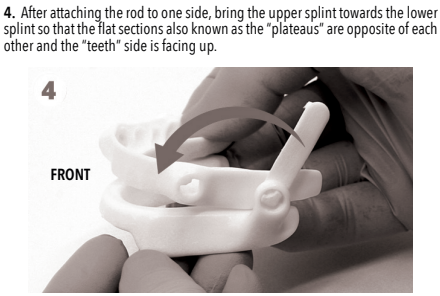Once both rods are in place, the plateaus (flat sides) of the splints should be facing each other, and the D-SAD is ready to use.

Advancing the D-SAD Classic for Optimal Treatment


During the advancement phase, the rods may be changed by the dentist to obtain the optimum therapy position of the D-SAD Classic. The dental provider will inform you if this is needed. To increase protrusion, use shorter rods (ex: go from 26 to 25.3 mm). To decrease protrusion, use longer rods (ex: go from 25.5 to 26 mm). Once the final setting is achieved, the rods are fixed in place and do not need to be changed. The anchors and the rods become the permanent fixed mechanical hinge of the appliance. Always use the same length of rods on both sides, unless specified by your dental practitioner.

Removing the rods during the advancement phase

Recreate a backward "L" with the device as noted above (Figure 5). Once in this position, turn the rod to release from the upper splint anchor. Then remove the other side.

D-SAD X3

- Parts
- Upper and lower splints.
 - A snap-in clip kit for jaw advancement.

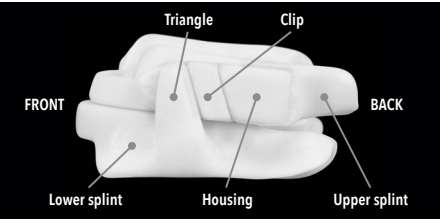
Note: The design of the device may vary depending on the specifics requested by your dental professional.

Clip Installation

- To insert the clip, push the open end of the clip in the clip housing opening.
- To remove the clip, hold the upper and lower part of the clip and pull out from the clip housing opening.

During the advancement phase, the clips may be changed by the dentist to obtain the optimum therapy position of the D-SAD X3. The dental provider will inform you if this is needed. To increase protrusion, use longer clips (ex: go from 0 to 0.5 mm). To decrease protrusion, use shorter clips (ex: go from 1 to 0.5 mm). Once the final setting is achieved, the clips are fixed in place and do not need to be changed. Always use the same length of clips on both sides, unless specified by your dental practitioner.

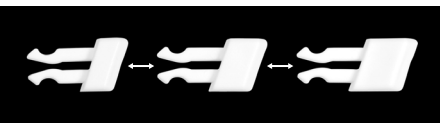Additional pictures are available at www.pantheradental.com/downloads

The next section applies to both devices

Getting ready for sleep

- Brush and floss your teeth and moisten your device with lukewarm water.
- As you place the device in your mouth, position the UPPER SPLINT FIRST on the upper teeth by pressing on both sides with your thumbs.
- Move your jaw forward, press the lower splint in place onto your lower teeth. Press on both sides with your fingers.
- Slowly close your mouth. In this configuration, the jaw is at a forward position.

Removing the D-SAD Classic or X3

- Remove the lower splint: place the thumbs on both sides of the lower splint, near the gum in front at the level of lower canines, and lift both sides of the device at the same time.
- Remove the upper splint: place your fingers on both sides of the upper splint, near the gum, and pull down on both sides of the splint at the same time.
- Remove the device from your mouth. Clean it according to the cleaning instructions.

Cleaning and storing instructions

- Every morning after each use, rinse your device in lukewarm water. Clean it with a soft bristle toothbrush or denture brush. Brush all surfaces and sides.
- Inspect your device for cleanliness and damage.
- Dry your device thoroughly with clean paper or cloth towel.
- Before storing your device, make sure the storage container is clean and dry.
- Clean your device at least twice a week with an antibacterial chlorine-free denture cleaner solution. Follow the manufacturer's instructions. Soak and brush the device with the solution. Rinse and dry well with a clean paper or cloth towel.

- Store in the storage container with the lid off, at room temperature away from direct heat, humidity and light. Never store your device in a closed container when it is humid.
- For extra cleaning, clean your device at least twice a week with an antibacterial chlorine-free denture cleaner solution. Follow the manufacturer's instructions. Soak and brush the device with the solution. Rinse and dry well with a clean paper or cloth towel.

Warnings

- Do not put your device in boiling water (100°C or 212°F).
 - Verify that you have no allergy to the ingredients present in the antibacterial chlorine-free denture cleaning solution before using it.
 - Do not use mouthwash to clean your device, this could stain and induce odor or taste.
 - Do not use alcohol-based products, bleach, acetic acid, household cleaning products or chlorine based solutions.
 - Note: Discoloration of the device may be a normal occurrence; this is due to the acidity of saliva. Some medicine, amalgam fillings, coffee, tea, cola, wine may affect the color of the device. This will not affect the performance or the warranty of the device.
- Helpful morning exercises
- If you are having any difficulty with your teeth in the morning:
 - Your dentist may provide instructions for morning exercise to help re-align your bite after wearing the device, such as moving the chin backward, clench molars together and hold for 2-3 seconds (repeat 5-10 times).
 - Eat breakfast or chew sugar-free gum; it will help your bite go back to normal.

Other tips

- Drink only water while the device is in place.
- Air Travel: It is advised to pack the device in your carry on.
- Keep away from pets as they may think it is a "chew toy".

Material

Both the D-SAD Classic and D-SAD X3 are made from 100% polyamide type 12, a medical grade material in the nylon family. Allergy and intolerance to this biocompatible material are extremely rare. If you have doubts, stop using the device and consult your dental professional or physician.

Warranty

The Panthera Sleep devices are guaranteed for 3 years from the date of manufacture. The warranty covers the devices against any workmanship or materials defects. This is a manufacturer's warranty and any repairs should only be done by Panthera Dental or their authorized partners/distributors. This warranty does not include any costs associated with additional dental or medical visits.

Note: To honor this warranty, the device and/or the rods or clips must be returned by your dental health care provider to Panthera Dental. Panthera Dental has the discretion to replace or repair the device.

USER ASSISTANCE INFORMATION

Your Dentist Phone number:

- For more information, visit www.pantherasleep.com
- Panthera Dental, Canada, Head Office
9105 rue John Simons, Quebec city, Quebec, Canada, G2B 0S6

Panthera Dental, France :
38 Rue Jean Mermoz, 78600 Maisons-Laffitte France.
Panthera Dental, USA :
22343 La Palma Ave, Ste 128, Yorba Linda, CA, USA 92887.

Validity:
The publication of these instructions for use cancels and replaces all previous versions. Authorization and availability of products may be limited to certain countries/regions. © Panthera Dental, 2022. All rights reserved.

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